Quick Reference Guide

Vivago

www.vivago.com

FIN

CFI9068-v1

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Vivago

www.vivago.fi

ITA

CIT9068-v1

Guida Rapida

Vivago

www.vivago.com

IF YOU THINK YOU NEED HELP, PLEASE PRESS THE WRISTCARE MANUAL ALARM BUTTON IMMEDIATELY TO RAISE AN ALARM. PRESS THE BUTTON UNTIL THE RED LIGHT OF THE WRIST UNIT GOES ON.

DO NOT WAIT FOR THE SYSTEM TO DO IT FOR YOU AUTOMATICALLY.

- The WristCare system sends an automatic alarm e.g. in case of immobility within 30 minutes to 4 hours. It is possible to adjust the sensibility of the automatic alarm to be user-specific.

- After an alarm has been raised, the alarm recipient will be able to talk with you via the base unit's loudspeaker (Voice communication is a separate function that is used in some of the systems)

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IMPORTANT!

WristCare can only raise an alarm when you are within the reach of the base station.

WristCare system can only send an alarm when the unit is properly fastened around your wrist.

WristCare cannot send an automatic alarm when you have taken the unit off your wrist or if there is fabric or other material between the wrist unit and your skin.

The WristCare system can only help you by sending an alarm, either manually or automatically. **DO NOT WAIT FOR THE SYSTEM TO RAISE AN ALARM FOR YOU AUTOMATICALLY IF YOU ARE ABLE TO DO IT YOURSELF.** Make sure you understand the purpose of this device and the limitation of the added safety it provides.

The WristCare system isn't a medical device. It shouldn't be used to monitor a specific medical condition.

IT IS IMPORTANT TO KEEP THE WRISTCARE ON YOUR WRIST AT ALL TIMES, EVEN AT NIGHT.

You should remove your WristCare unit from your wrist in order to:

- Go to the sauna or other hot environments.
- Clean the WristCare unit once a week according to the instructions in the User's Guide.

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**The button functions**

WristCare has three buttons:

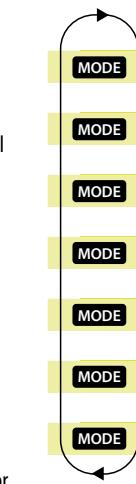
1 VIVAGO button on the watch is meant for triggering the alarms, except when doing confirming changes in the setting mode of the watch.

2 A long press on LIGHT button will switch on the background light. The light will automatically switch off after 12 seconds.

3 MODE button is to be used when browsing through the different menus.

By pressing LIGHT(2) and MODE(3) buttons in the watch mode at the same time you enter the set-up mode, where you can set the time and date

A long press of the MODE button or no press for 60 seconds will, from all menus, redirect the system to the watch mode.

The Menu Structure

Settings
By pressing and holding both the LIGHT and MODE buttons simultaneously you enter the setting mode. The setting mode depends on in which display mode you enter the settings.

In the set up mode the buttons work as follows:
LIGHT: changes settings.
MODE: changes the numbers, for instance from hours to minutes when setting the time or changes the settings.

VIVAGO: confirms the setting and moves you to next set-up mode.

If you don't press any button for 60 seconds you will automatically return to the basic watch mode

Please notice the correct positioning of the charger.

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